



ALL IN ONE ELECTRICAL MOBILE HOIST

For comfort and safety

TECHNICAL SPECIFICATIONS, PATIENT HOIST

		" See page 4
Safe working load	150 kg/330 lb.	* 200 kg/440 lb.
Spreader bar - High Position	180 cm/70.86 in.	177 cm/69.68 in.
Spreader bar - Low position	55 cm/21.65 in.	58 cm/22.84 in.
Lower chassis frame, length	117 cm/46.06 in.	122 cm/48.03 in.
Overall height of legs	12 cm/4. 72 in.	12 cm/4.72 in.
Width - legs closed	65 cm/25.59 in.	69 cm/27.16 in.
Width - legs opened	109 cm/42.91 in.	133 cm/52.36 in.
Ball bearing front castor	7,5 cm/2.95 in.	7,5 cm/2.95 in.
Rear castors with brake	10 cm/3.93 in.	10 cm/3.93 in.
Battery Box weight	3 kg/6.6 lb.	3 kg/6.6 lb.
Total weight	51 kg/112 lb.	54,3 kg/119.5 lb.
Approx. number of lifts per charge 60		30



Recharging is simply done by plugging the power cord into the wall socket.



The unique "All-in-One" active Hoist can be converted from a standard Patient Sling Hoist to a Stand-Up Hoist - without using any tools. It is easy to convert the Patient Hoist to a special Stand-Up Hoist and make transfers in a standing position.



With the "All-in-One" hoist it is possible to lift a client from the floor.



The Patient Hoist is simple to operate and has a wide lifting area - from lowest to highest position - facilitating floor to high bed lifts.

- The "All-in-One" hoist has a very small turning radius, as the chassis incl. castors has a length of only 117 cm/46.06 in.
- The "All-in-One" hoist is compact and practical but still combines strength with simplicity in use.
- The "All-in-One" hoist is safe and easy to operate, and incorporates many unique safety features, e.g., the hoist will stop if the lifting bar hits an obstruction.
- Notice, that the height of the lower chassis is only 12 cm/4.72 in. which means that the hoist can easily be positioned under most beds.



The special Domino "Hygiene" sling provides excellent access to the client, it gives good support and is easy to use.



TECHNICAL SPECIFICATIONS, STAND-UP HOIST

		^ Se side 4
Safe working load	150 kg/330 lb.	* 200 kg/440 lb.
Fork (Bar) - High position	191 cm/75.20 in.	191 cm/75,2 in.
Fork (Bar) - Low position	95 cm/37.40 in.	95 cm/37.40 in.
Lower chassis frame, length	117 cm/46.06 in.	122 cm/48.03 in.
Overall height of legs	12 cm/4.72 in.	12 cm/4.72 in.
Width - legs closed	65 cm/25.59 in.	69 cm/27.16 in.
Width - legs opened	109 cm/42.91 in.	133 cm/52.36 in.
Ball bearing from castors	7,5 cm/2.95 in.	7,5 cm/2.95 in.
Rear castors with brake	10 cm/3.93 in.	10 cm/3.93 in.
Battery box weight	3 kg/6.6 lb.	3 kg/6.6 lb.
Standing platform &		
knee support weight	8,5 kg/18.7 lb.	8,5 kg/18.7 lb.
Total weight	58,5 kg/128.7 lb.	61,5 kg/135.3 lb.
Approx. numbers of lifts per	charge 60	30

When the Patient Hoist is converted into a Stand-Up Hoist, the result is an Active Hoist. The client now becomes an active part of the entire lifting process, during which the client is transferred in a dignified manner. The client experiences an advantage in being transported in a standing position.



Notice: The knee support is anatomically designed - and can be adjusted both vertically and horizontally. That means you always achieve a correct adjustment of the knee support.

When a change of the Patient Hoist to a Stand-Up Hoist is required, the outer part of the lifting bar is removed and the fork shaped bar + foot and knee support is attached with no use of tools, at all. For the sake of convenience the hoist should be braked during conversion.



Effortless raising and lowering of client. Easy and comfortable transfer between wheelchairs, beds, toilet, bath, etc. Clear access to the clients clothing for toileting. Elimination of stress and strain upon the carers back. The pictures below shows how easy a client can be transferred from a wheelchair to a toilet.



Being lifted is pleasant for the client. No pressure is experienced under the arms.



The client can gently be lowered on to the toilet.



The hoist is equipped with both electric and manual operated emergency lowering. The battery is easily interchangeable.


ORDERING INFORMATION

ITEM NO. DESCRIPTION

25.20020	COMBINATION HOIST ALL-IN-ONE 150 KG/330 LB.			
25-20020	All-in-One combined Patient & Stand-Up Hoist. Can be altered from a			
	Patient Hoist to a Stand-up Hoist and vice versa. Incl. battery and charger.			
	PATIENT HOIST 150 KG/330 LB.			
25-20040	Patient Hoist, standard, incl battery and charger. *			
25-20010	Patient Hoist, standard, incl battery and charger.			
	Can be altered to a Stand-up Hoist by use of Item. No. 25-20120			
	STAND-UP HOIST 150 KG/330 LB.			
25-20050	Stand-up Hoist, standard, incl battery and charger. *			
25-20030	Stand-up Hoist, standard, incl battery and charger.			
	Can be altered to a Patient Hoist by use of Item. No. 25-20455			
	COMBINATION HOIST ALL-IN-ONE 200 KG/440 LB.			
25-20520	All-in-One combined Patient & Stand-Up Hoist. Can be altered from a			
	Patient Hoist to a Stand-up Hoist and vice versa. Incl. battery and charger.			
	PATIENT HOIST 200 KG/440 LB.			
25-20524	Patient Hoist, standard, incl battery and charger. *			
25-20521	Patient Hoist, standard, incl battery and charger.			
	Can be altered to a Stand-up Hoist by use of Item. No. 25-20517			
	STAND-UP HOIST 200 KG/440 LB.			
25-20525	Stand-up Hoist, standard, incl battery and charger. *			
25-20523	Stand-up Hoist, standard, incl battery and charger.			
	Can be altered to a Patient Hoist by use of Item. No. 25-20516			
	ACESSORIES			
25-20120	Stand-up kit, All-in-One 150 kg/330 lb. Fork shaped lifting bar, standing			
	platform with adjustable knee support. For Item No. 25-20010			
25-20517	Stand-up kit, All-in-One 200 kg/440 lb. Fork shaped lifting bar, standing			
	platform with adjustable knee support. For Item No. 25-20521			
25-20455	Patient Hoist kit, All-in-One 150 kg/330 lb. Lifting bar and spreader bar.			
	For Item No. 25-20030			
25-20516	Patient Hoist kit, All-in-One 200 kg/440 lb. Lifting bar and spreader bar.			
	For Item No. 25-20523			
25-97260	Padding for spreader bar.			
25-29950	Wall charging depot. Incl. battery.			
25-29030	Battery.			
25-29235	Digital scale. Capacity 300 kg/660 lb.			
	* Is fitted with fixed arm			

* Is fitted with fixed arm.

Ask for the complete sling brochure with detailed description of Domino Slings.

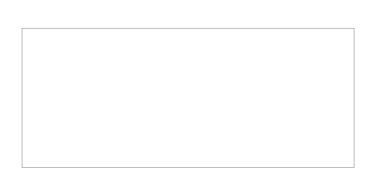







Ropox A/S Ringstedgade 221 4700 Naestved Telefon 5575 0550 Telefax 5575 0550 info@ropox.dk www.ropox.COM

1962349 0321.